

## PATIENT INFORMATION

### Follistim® (Fol'-lis-tim) AQ Cartridge (follitropin beta injection)

Read the Patient Information that comes with Follistim® AQ Cartridge before you start using it and each time you get a refill. There may be new information. This information does not take the place of talking with your healthcare provider about your medical condition or treatment.

#### What is Follistim AQ Cartridge?

Follistim AQ is a prescription medicine that contains follicle-stimulating hormone (FSH). The medicine is taken with the Follistim Pen®.

#### Follistim AQ Cartridge is used:

##### In women:

- to help healthy ovaries to develop (mature) and release eggs
- as part of treatment programs that use special techniques (skills) to help women get pregnant by causing their ovaries to produce more mature eggs

##### In men:

- to help bring about the production and development of sperm

#### Who should not take Follistim AQ Cartridge?

##### Do not take Follistim AQ Cartridge if you are a Woman or Man who:

- is allergic to recombinant human FSH products
- has a high level of FSH in your blood indicating that your ovaries (women only) or testes (men only) may be permanently damaged and do not work at all
- has uncontrolled thyroid, pituitary, or adrenal gland problems
- is allergic to streptomycin or neomycin (types of antibiotics)
- has a tumor of the hypothalamus, pituitary gland, breast, uterus (women only), ovary (women only), or testis (men only)

##### Do not take Follistim AQ Cartridge if you are a Woman who:

- is pregnant or think you may be pregnant
- has heavy or irregular vaginal bleeding and the cause is not known
- has ovarian cysts or enlarged ovaries, not due to polycystic ovary syndrome (PCOS)

Talk to your healthcare provider before taking this medicine if you have any of the conditions listed above.

#### What should I tell my healthcare provider before taking Follistim AQ Cartridge?

##### Before you take Follistim AQ, tell your healthcare provider if you:

- have ever had Ovarian Hyperstimulation Syndrome (OHSS).
- are pregnant or think you may be pregnant. Follistim AQ Cartridge should not be used in pregnant women.
- have an increased risk of blood clots (thrombosis).
- have ever had a blood clot (thrombosis), or anyone in your immediate family has ever had a blood clot (thrombosis).
- have ever had stomach (abdominal) surgery.
- have ever had twisting of your ovary (ovarian torsion).
- have or have had a cyst in your ovary.
- have polycystic ovarian syndrome (PCOS).
- have any other medical conditions (for example, diabetes, heart disease, or any other long-term disease).
- are breastfeeding or plan to breastfeed. It is not known if the medicine in Follistim AQ Cartridge passes into your breast milk. You and your healthcare provider should decide if you will take Follistim AQ Cartridge or breastfeed. You should not do both at the same time.
- have been told by a doctor that pregnancy would be dangerous for you.

**Tell your healthcare provider about all the medicines you take**, including prescription and non-prescription medicines, vitamins, and herbal supplements.

Know the medicines you take. Keep a list of them and show your healthcare provider and pharmacist when you get a new medicine.

#### How should I use Follistim AQ Cartridge?

- Be sure that you read, understand, and follow the "Patient Instructions for Use" that come with Follistim AQ Cartridge.
- Use Follistim AQ Cartridge exactly as your healthcare provider tells you to.
- Your healthcare provider will tell you how much Follistim AQ Cartridge to use, how to inject it, and how often it should be injected.
- Do not inject Follistim AQ Cartridge at home until your healthcare provider has taught you the right way to put the cartridge and pen device together and to inject yourself.
- Do not mix any other medicines into the cartridge.
- Do not change your dose of Follistim AQ Cartridge unless your healthcare provider tells you to.
- Call your healthcare provider immediately if you use too much Follistim AQ Cartridge.
- If you miss or forget to take a dose, do not double your next dose. Ask your healthcare provider for instructions.
- Use Follistim AQ Cartridge only with the Follistim Pen.
- Do not use the Follistim Pen if you are blind or visually impaired unless you have assistance from an individual with good vision who is trained in the right way to use the pen.

- Do not re-use the BD Micro-Fine™ Pen Needle.
- Your healthcare provider may do blood hormone tests while you are taking Follistim AQ Cartridge. Make sure you follow-up with your healthcare provider to have your blood tested when told to do so.

**Women:**

- Your healthcare provider will do ultrasound scans of your ovaries to monitor how the follicles are growing. Make sure you follow-up with your healthcare provider to have your ultrasound scans.

**Men:**

- Your healthcare provider may test your semen while you are taking Follistim AQ Cartridge. Make sure you follow-up with your healthcare provider to give a semen sample for testing.

**What are the possible side effects of Follistim AQ Cartridge?**

**Follistim AQ Cartridge may cause serious side effects.**

**Serious side effects in women include:**

- Ovarian enlargement. Your ovaries may become very large at the end of your fertility treatment cycle. Therefore, close supervision by your healthcare provider is very important. To check the effects of your fertility treatment, your healthcare provider will take ultrasound scans of your ovaries and may also check blood hormone levels.
- Ovarian hyperstimulation syndrome (OHSS). OHSS is a serious medical problem that can happen when the ovaries are over stimulated. In rare cases, OHSS may be life-threatening and has caused death. OHSS causes fluid to build up suddenly in your stomach and chest areas and can cause blood clots to form. Call your healthcare provider right away if you have:
  - pain in your lower stomach area.
  - nausea.
  - vomiting.
  - sudden weight gain due to fluid buildup.
  - diarrhea.
  - decreased urine output.
  - trouble breathing.
- Lung problems. Treatment with gonadotropic hormones like Follistim AQ Cartridge can cause you to have trouble breathing (acute respiratory distress syndrome).
- Blood clots. Treatment with gonadotropic hormones like Follistim AQ Cartridge may (just as with pregnancy) increase your chance of having blood clots (thrombosis) in your blood vessels. Thrombosis is the formation of a blood clot in a blood vessel. Blood clots can lead to:
  - blockage in your lungs (pulmonary embolus).
  - stroke.
  - heart attack.
  - blood vessel problems (thrombophlebitis).
  - a lack of blood flow (deep venous thrombosis) that may result in a loss of your arm or leg.
- Ovarian torsion. Follistim AQ Cartridge may increase the chance of twisting of the ovaries in women with certain conditions such as OHSS, pregnancy and previous stomach (abdominal) surgery, past history of a twisting of an ovary, past or current cysts in your ovary or ovaries, or polycystic ovaries. Twisting of the ovary could cause the blood flow to the ovary to be cut off.
- Pregnancy and birth of multiple babies. There is an increased chance of having twins or more than two babies, even when only one embryo is transferred into the uterus. Having a pregnancy with more than one baby at a time increases the health risk for you and your babies. Discuss your chances of multiple births with your healthcare provider.
- Birth defects. A woman's age, certain sperm problems, genetic background of both parents and a pregnancy with multiple babies can increase the chance that your baby might have a birth defect.
- Ectopic pregnancy (pregnancy outside of the uterus). If treatment with Follistim AQ Cartridge results in pregnancy, there is a higher chance of pregnancy outside the uterus. Therefore, your healthcare provider should perform an early ultrasound examination to exclude the possibility of pregnancy occurring outside of the uterus.
- Miscarriage. Miscarriages have been reported for all gonadotropin treatments, including Follistim AQ Cartridge. However, it is not known if miscarriages are related to treatment with fertility medicines or the underlying infertility.
- Ovarian and Other Reproductive System Tumors. There have been reports of ovarian and other reproductive system tumors in women who have had infertility treatment. It is not known if treatment with fertility medicines increases the risk of these tumors in infertile women.

**The most common side effects of Follistim AQ Cartridge include:**

**In women:**

- headache
- nausea
- stomach pain
- discomfort or pain in the lower stomach area
- cyst (closed sac) in the ovary
- feeling tired

**In men:**

- headache
- pain at the injection site
- bruising, swelling or redness at the injection site
- breast enlargement
- acne

These are not all the possible side effects of Follistim AQ Cartridge. For more information, ask your healthcare provider or pharmacist.

**Call your healthcare provider immediately** if you get worsening or strong pain in the lower stomach area (abdomen). Also, call your healthcare provider immediately if this happens some days after the last injection has been given.

Tell your healthcare provider if you have any side effect that bothers you or that does not go away.

Call your doctor for medical advice about side effects. You may report side effects to FDA at 1-800-FDA-1088.

#### **How should I store Follistim AQ Cartridge?**

- Store Follistim AQ Cartridge in the refrigerator between 2°-8°C (36°-46°F) until the expiration date.
- Follistim AQ can be stored at or below 25°C (77°F) for 3 months or until the expiration date, whichever comes first. Once the rubber inlay of the Follistim AQ Cartridge has been pierced by a needle, the product may be stored only for a maximum of 28 days at 2°-25°C (36°-77°F).
- Keep Follistim AQ Cartridge away from light.
- Do not freeze.

**Keep Follistim AQ Cartridge and all medicines out of the reach of children.**

#### **General information about Follistim AQ Cartridge**

Medicines are sometimes prescribed for purposes other than those listed in the Patient Information leaflet. Do not use Follistim AQ for a condition for which it was not prescribed. Do not give Follistim AQ Cartridge to other people, even if they have the same condition that you have. It may harm them.

This Patient Information leaflet summarizes the most important information about Follistim AQ Cartridge. If you would like more information, talk with your healthcare provider. You can ask your pharmacist or healthcare provider for more information about Follistim AQ Cartridge that is written for healthcare professionals.

For more information, go to [www.follistim.com](http://www.follistim.com) or call 1-866-836-5633.

#### **What are the ingredients in Follistim AQ Cartridge?**

Active ingredient: follitropin beta

Inactive ingredients: sucrose, sodium citrate, benzyl alcohol NF-10 mg/mL, L-methionine, polysorbate 20, water for injection, hydrochloric acid, and/or sodium hydroxide.

Manufactured for: Merck Sharp & Dohme Corp., a subsidiary of  
 **MERCK & CO., INC.**, Whitehouse Station, NJ 08889, USA

Manufactured by: Vetter Pharma-Fertigung GmbH & Co. KG, Ravensburg, Germany

BD, BD Logo and BD Micro-Fine are trademarks of Becton, Dickinson and Company.

For patent information: [www.merck.com/product/patent/home.html](http://www.merck.com/product/patent/home.html)

Copyright © 2004-2019 Merck Sharp & Dohme B.V., a subsidiary of **Merck & Co., Inc.**  
All rights reserved.

Revised: 09/2019

usppi-mk8328-SOI-1909r014

## PATIENT INSTRUCTIONS FOR USE

### Follistim® (Fol'-lis-tim) AQ Cartridge (follitropin beta injection)

Read the Patient Instructions for Use that comes with Follistim® AQ Cartridge before you start using it and each time you get a refill. There may be new information. This information does not take the place of talking with your healthcare provider about your medical condition or treatment.

#### A. Getting Ready

- Follistim Pen is not recommended for the blind or visually impaired user without the assistance of an individual with good vision, trained in the proper use of the injection device.
- Learn about all of the parts of the Follistim Pen (See Figure 1), Follistim AQ Cartridge (See Figure 2) and the BD Micro-Fine™ Pen Needle (See Figure 3). You will need to recognize these parts to follow the directions.



Figure 1. Follistim Pen and its Parts

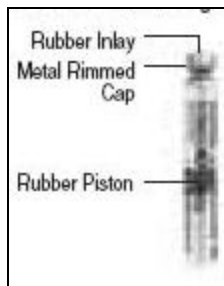


Figure 2. Parts of Follistim AQ Cartridge

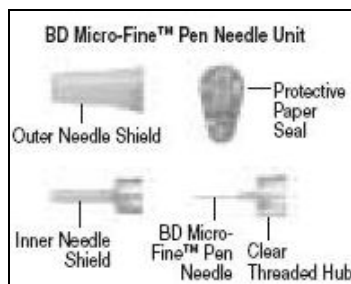


Figure 3. Parts of BD Micro-Fine Pen Needle Unit

- Remove the Cartridge out of the refrigerator.
- Injecting cold drug is likely to cause discomfort. Therefore, it is recommended you allow the drug to reach room temperature before taking the injection.
- Check the liquid in the cartridge. It should appear clear and colorless. If the solution is not clear and colorless or has particles in it, **do not use it**.
- **Gather the supplies you will need for your injection. You will need:**
  - a clean dry surface
  - alcohol
  - cotton balls or alcohol pads
  - sterile gauze
  - a puncture-proof container to throw away the used syringe and needle
- Wash your hands with soap and water and dry them before you use Follistim Pen or when you replace the cartridge.


**B. Loading the Follistim Pen with the Follistim AQ Cartridge**

- Holding the Pen Body firmly with one hand, pull off the Protective Cap with your other hand (See Figure 4). Put the cap aside on a clean, dry surface.




**Figure 4**

- Unscrew the entire Pen Body from the Cartridge Holder (See Figure 5). Place the Cartridge Holder and the Pen Body aside on a clean, dry surface.

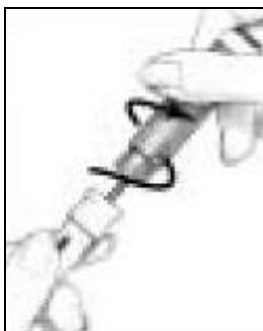

**Figure 5**

- Take the Follistim AQ Cartridge out of its package. Clean the rubber inlay on the cartridge with an alcohol pad. Pick up the Cartridge Holder. Put the Cartridge into the Cartridge Holder (See Figure 6). The metal rimmed cap goes in first.

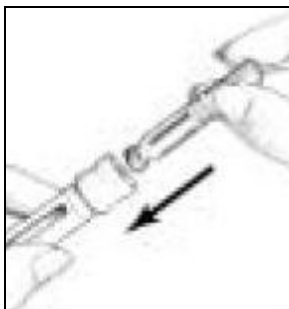

**Figure 6**

- Pick up the Pen Body and lower it into the Cartridge Holder. The black rod must press against the Rubber Piston on the cartridge. Screw the Pen Body fully onto the Cartridge Holder (See Figure 7). Make sure there is no gap between the Pen Body and the Cartridge Holder. The arrow (▲) on the Cartridge Holder should point to the middle of the yellow alignment mark (■) on the blue Pen Body.



**Figure 7**

- Clean the open end of the Cartridge Holder with an alcohol pad (See Figure 8).

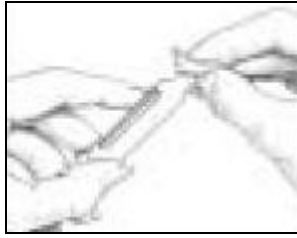


Figure 8

- Pick up a new BD Micro-Fine Pen Needle that is in its Outer Needle Shield. Peel off the protective paper seal (See Figure 9). Do **not** touch the needle. Do not place the open needle on any surface. **Use Only the BD Micro-Fine 0.33 mm × 12.7 mm (29G) Pen Needles as supplied with the Follistim AQ Cartridge.**
- You must use a new BD Micro-Fine Pen Needle with each injection. Never reuse a needle. Attach a new BD Micro-Fine Pen Needle after you make sure there is a Follistim AQ Cartridge in the Cartridge.

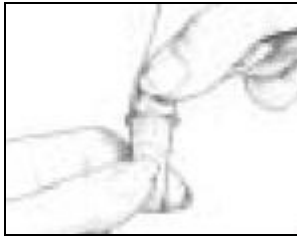


Figure 9

- Hold the Outer Needle Shield firmly in one hand while holding the Cartridge Holder firmly in the other hand. Push the end of the Cartridge Holder into the Outer Needle Shield. Screw them tightly together (See Figure 10). Place your Follistim Pen with the loaded cartridge and attached needle, flat on a clean, dry surface.

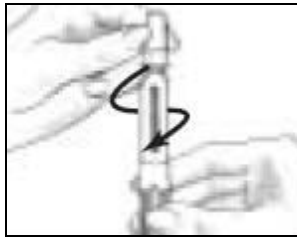


Figure 10

### C. Preparing the Injection Site

- Follistim AQ Cartridge can be injected directly into a layer of fat under your skin (subcutaneously).
- When giving a subcutaneous injection, follow your healthcare provider's instructions about changing the site for each injection. This will help lower your chances of having a skin reaction.
- **Do not** inject Follistim AQ Cartridge into an area that is tender, red, bruised, or hard.
- The recommended site for injecting Follistim AQ Cartridge subcutaneously is:
  - just below your belly button (navel) (See Figure 11)



Figure 11

- the upper outer area of your thigh (See Figure 12)

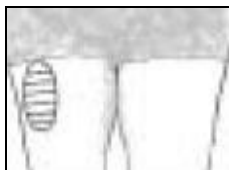


Figure 12

- Clean the skin with an alcohol wipe where the injection is to be made. Clean about two inches around the injection site where the needle will be inserted. Do not touch the cleaned area of skin.

#### D. Dialing the Dose Before You Give the Injection

- Your healthcare provider will decide on the dose of Follistim AQ Cartridge to be given. This dose may be increased or decreased as your treatment progresses depending on your individual type of treatment.
- Follistim AQ Cartridge using Follistim Pen can be administered subcutaneously (beneath the skin) in prescribed doses from 50 International Units (IU) up to 450 IU, in marked 25 IU increments. The Dosage Scale on the Pen has numbers and audible clicks to help you set the correct dose.
- Pull off the outer needle shield. Leave the Inner Needle Shield in place over the needle attached to the Pen (See Figure 13). Do not throw the Outer Needle Shield away, you will need it later when you throw the needle away.

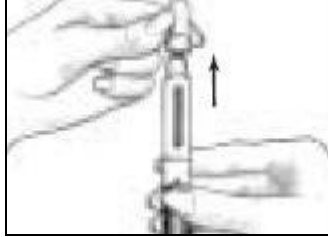


Figure 13

- Carefully remove the Inner Needle Shield and discard it (See Figure 14). Do not touch the needle or let it touch any surface while uncapped.



Figure 14

- Hold the Follistim Pen with the needle pointing upwards. Tap the Cartridge Holder gently with your finger to help air bubbles rise to the top of the needle. The small amount of air bubble will not affect the amount of medicine you receive.
- With a loaded new, unused cartridge:
  1. Dial the Dosage Knob until you hear one click. With the needle pointing upwards, push in the Injection Button.
  2. Look for a droplet at the tip of the needle (See Figure 15). If you see the droplet, then you can dial in your dose.
  3. If you do **not** see a droplet, repeat Step 1 (as above) until you see droplet.

You must **make sure you see a droplet** of medicine (**check the flow of medicine**) or you may **not** inject the correct amount of medicine.



Figure 15

- With a partially used cartridge, to give yourself another dose of medicine you will need to attach a new BD Micro-Fine Pen Needle and look for a droplet forming at the tip of the needle (See Figure 15 above). If you see a droplet, then you can dial in your dose.
- If **no** droplet:
  1. Dial the Dosage Knob until you hear one click. With the needle pointing upwards, push in the Injection Button.
  2. Look for a droplet at the tip of the needle. If you see the droplet, then you can dial in your dose.
- Your Follistim AQ Cartridge should be one of the following:
  - Orange Metal Cap – 150 international units
  - Silver Metal Cap – 300 international units

- Gold Metal Cap – 600 international units
- Blue Metal Cap – 900 international units

If you did **not** understand that you should have one of the cartridges above, please contact your healthcare provider.

- For doses of 50 IU up to 450 IU, turn the Dosage Knob until the correct dosage aligns with the dosage markers on each side of the Dosage Window (see Figure 16).



Figure 16

- If by mistake you dial past the correct number, do not try to turn the Dosage Knob backward to fix the mistake. Continue to turn the Dosage Knob in the same direction past the 450 IU mark, as far as it will turn. The Dosage Scale must move freely. Push the Injection Button in all the way. See Figure 17. Start to dial again starting from "0" upwards. By following these directions, you will not lose any medicine from the Follistim AQ Cartridge (See "Checking the Medicine Level Remaining").

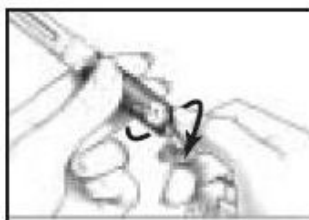


Figure 17

- If you turn the Dosage Knob backward to correct the mistake, it will not damage the Pen, but you will lose some medication from the Follistim AQ Cartridge.
- Never dial your dose or try to correct a dialing mistake when the needle is still in your skin as this may result in your receiving an incorrect dose.
- If your prescribed dose exceeds the deliverable dose of Follistim Pen or exceeds the amount remaining in the cartridge, you will need to give yourself more than one injection.

#### E. Giving Yourself an Injection

- Pinch a fold of skin at the cleaned injection site. **Do not** touch the cleaned area of skin.
- With the other hand hold the entire Pen with Cartridge loaded and Needle on like you would a pencil. Use a quick "dart-like" motion to insert the needle straight up and down (90-degree angle).
- Press the injection button all the way in to make sure you give yourself a full injection. (See Figure 18). Wait for five seconds before pulling the needle out of the skin. The middle of the Dosage Window should display a dot next to the "0".





**Figure 18**


If the injection button does **not** push in all the way, and the number in the Dosage Window does not read "0", it means there is not enough medication left in the cartridge to complete your prescribed dose. The number in the Dosage Window will give you the amount of medicine needed to complete your dose. Write this number down. This will be the number you dial for the completion of your dose. **Start over** with a new Follistim AQ Cartridge and a new needle and follow all the instructions up to this step. Make sure you choose a different injection site to complete your dose of Follistim AQ Cartridge.

- Pull out the BD Micro-Fine Needle and firmly press down on the injection site with an alcohol swab. Use the BD Micro-Fine Pen Needle for one injection only.
- Place the Outer Needle Shield on a flat table surface with the opening pointing up. The opening of the Outer Needle Shield is the wider end with the rim. Without holding on to the Outer Needle Shield, carefully insert the needle (attached to the Follistim Pen) into the opening of the Outer Needle Shield and push down firmly. The Outer Needle Shield should now be attached to the Cartridge Holder and covering the needle (See Figure 19).



**Figure 19**

- Grip the Outer Needle Shield and use it to unscrew the needle from the Cartridge Holder (See Figure 20). If there is Follistim AQ Cartridge medicine left for more injections, put the Pen Cap back on the Pen Body and store your Follistim Pen in a safe place in the refrigerator (not in the freezer) or at room temperature. Never store the Follistim Pen with a needle attached to it. If you are giving an injection to another person, be very careful when removing the needle from the skin. Accidental needle sticks can transmit potentially serious or grave infectious diseases.

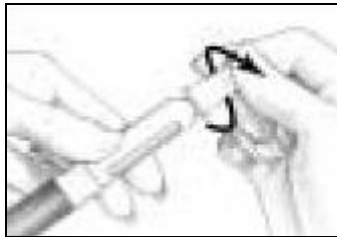

**Figure 20**

- Throw away the Outer Needle Shield with the used needle right away. Do not throw away in a trash can. Place it in a "special" container. (See "How Do I Throw Away Used Cartridges and Needles?")
- If there is Follistim AQ Cartridge medicine left for more injections, put the Pen Cap back on the Pen Body and store your Follistim Pen in a safe place in the refrigerator (not in the freezer) or at room temperature. Never store the Follistim Pen with a needle attached to it. If you are giving an injection to another person, be very careful when removing the needle from the skin. Accidental needle sticks can lead to serious infections.
- Unscrew the Pen Body from the Cartridge Holder with the **empty** Follistim AQ Cartridge (See Figure 21).



**Figure 21**

- Put the Pen Body down on a clean, dry surface and remove the empty Follistim AQ Cartridge from the Cartridge Holder (See Figure 22). Safely, dispose of the empty Follistim AQ Cartridge right away in the same "special" container that you used for the needle disposal. Do not put the cartridge in a trash can. At the end of your treatment cycle, your doctor can advise you on how to properly dispose of the container. (See "How Do I Throw Away Used Cartridges and Needles?")

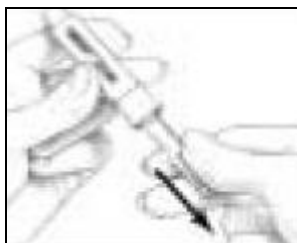


Figure 22

## F. Checking the Medicine Level Remaining

### For women and men:

Your healthcare provider should advise you of the number of prescribed doses which can be extracted from the full unused Follistim AQ Cartridge.

- **Do not** use the cartridge beyond the advised number of doses. Otherwise, you will run the risk that there will not be enough volume of drug for your prescribed dose.

### For women only:

- Keep a Follistim Pen Treatment Diary as follows:
  1. Record the Follistim AQ Cartridge content on Day 1. This will either be 150, 300, 600 or 900 international units depending on what your healthcare provider has prescribed for you.
  2. Record the dose you have been prescribed for your injection.
  3. Subtract your Day 1 dose from the Follistim AQ Cartridge content (150, 300, 600 or 900 international units). (See example – Figure 23.) This will give you the **remaining** Follistim AQ Cartridge content after the Day 1 dose is taken.
  4. Place the number determined as the content after Day 1 (see number 3) in the box as the Follistim AQ Cartridge content **available** for Day 2. (See example – Figure 23.)
  5. Subtract your Day 2 dose from the Follistim AQ Cartridge content you recorded in Step 4. This will give you the **remaining** Follistim AQ Cartridge content after Day 2. Record this number of units. (See example – Figure 23.)
  6. Repeat the steps to determine the Follistim AQ Cartridge content **available** and Follistim AQ Cartridge **remaining** for each day of use.

| Day | Date           | Dose Prescribed | Follistim AQ Cartridge Content Available | Follistim AQ Cartridge Content Remaining |
|-----|----------------|-----------------|--|--|
| 1   | month/day/year | 150             | 600                                      | 450                                      |
| 2   | month/day/year | 150             | 450                                      | 300                                      |
| 3   | month/day/year | 150             | 300                                      | 150                                      |
| 4   | month/day/year | 150             | 150                                      | 0  |

Figure 23. Example of Treatment Diary Starting with a 600 International Unit Cartridge

If you do not know if there is not enough medicine left in the Follistim AQ Cartridge for your next prescribed dose, see section “If There is Not Enough Follistim AQ in the Cartridge”.

## G. If There is Not Enough Follistim AQ in the Cartridge

1. If you realize **before** you inject that you do not have enough medicine remaining in your Follistim AQ Cartridge for your complete dose, follow either Option 1 or Option 2, but **not** both.
  - Option 1:
    - Dial your dose and inject the remaining content in the Follistim AQ Cartridge. The Dosage Knob Injection Button will not push in all the way (do not try to force down the button) and the Dosage Window number will not read “0” but will read the number of units you will need to complete your prescribed dose.
    - Write down the number of units needed to complete your dose.
    - Remove the needle and dispose of it properly (see “How Do I Throw Away Used Cartridges and Needles?”).
    - Using the Dosage Knob, reset the Dial Window to “0” by turning the Dosage Knob past the 450 IU mark as far as it will turn and push the Injection Button in all the way.
    - Before attempting to replace a Follistim AQ Cartridge, be sure that a BD Micro-Fine Pen Needle is not attached to the Follistim Pen.
    - Insert a new cartridge into the Follistim Pen and attach a new BD Micro-Fine needle.
    - Dial to the number of units you have written down to complete your prescribed dose.
    - Prepare a different injection site and inject the remaining drug to complete your dose (See “Giving Yourself an Injection”).
  - Option 2
    - Remove the Follistim AQ Cartridge.
    - **Start over** with a new Follistim AQ Cartridge and Insert into the Follistim Pen.
    - Follow the instructions for “Dialing the Dose” and “Giving Yourself an Injection.”
2. If you realize **after** you have inserted the needle at the injection site that you do not have enough medicine remaining in your Follistim AQ Cartridge for your complete dose:
  - Inject the remaining content in the Follistim AQ Cartridge. The Injection Button will not push in all the way and the number in the Dosage Window will not read “0” but will read the number of units you will need to complete your prescribed dose.
  - Wait 5 seconds before withdrawing the needle from your skin and gently apply pressure to the injection site with an alcohol pad.
  - Dispose of the used needle (See “How Do I Throw Away Used Cartridges and Needles?”).
  - Write down the number of units needed to complete your dose.
  - Using the Dosage Knob, reset the Dial Window to “0” by turning the Dosage Knob past the 450 IU mark as far as it will turn and push the Injection Button in all the way.

- Insert a new cartridge into the Follistim Pen and attach a new BD Micro-Fine needle.
- Dial to the number you have recorded to complete your prescribed dose.
- Prepare a different injection site and inject the remaining drug to complete your dose (See “Giving Yourself an Injection”).

#### H. How to Solve Problems with Follistim AQ Cartridge and Follistim Pen

If you have problems with using the Follistim AQ Cartridge and the Follistim Pen, see the following chart. If you still have problems after following the chart or if your problem is not on the chart, contact your healthcare provider.

| PROBLEM  | POSSIBLE CAUSES   | WHAT TO DO  |
|--|---|---|
| The Pen Body will not screw tightly into the Cartridge Holder.   | Is something in the way?  | Take out the Follistim Cartridge and check the Cartridge Holder to see if anything is in the way. Follow the instructions in this pamphlet to Screw the Pen Body fully onto the Cartridge Holder.   |
| No drug is coming out while checking the flow.   | The Cartridge Holder and the Pen Body are not properly screwed together.                  | Remove the current needle; tighten the Pen Body to the Cartridge Holder ensuring the arrow on the Cartridge Holder is pointing to the middle of the yellow alignment mark on the blue Pen Body. Attach a new needle to the Pen.<br><br>Recheck the flow as follows:<br>a. Dial the Dosage Knob until you hear one click. With needle pointing upwards, push in the Injection Button.<br>b. Look for a droplet at the tip of the needle. |
|  | Is the Follistim Cartridge empty?   | Change to a new cartridge.  |
|  | Has the needle been properly attached to the Follistim Pen?                               | Remove needle and replace with a new one, ensuring that the needle is screwed on tightly to the Pen.<br><br>Recheck the flow as follows:<br>a. Dial the Dosage Knob until you hear one click. With needle pointing upwards, push in the Injection Button.<br>b. Look for a droplet at the tip of the needle.  |
| You are concerned that you can turn the Dosage Knob to the next number without clicking and the injection button spins freely.                     | This is not a problem.  | The system is in the reset mode. The Injection Button and Dosage Knob must be pushed all the way down to '0' to re-engage the mechanism and the correct dose can now be set. A click will be heard for each setting in the viewing window.  |
| The Dosage Knob does not go back to '0' while you are injecting.   | Is the Follistim Cartridge empty?   | Change to a new cartridge.  |
|  | Is the needle blocked?  | a. Remove the needle from the skin and dispose of safely.<br>b. Check the Dosage Window and note how much remaining drug to inject.<br>c. Attach a new needle.<br><br>Recheck the flow as follows:<br>a. Dial the Dosage Knob until you hear one click. With needle pointing upwards, push in the Injection Button.<br>b. Look for a droplet at the tip of the needle.<br>c. Dial remaining dose.                                       |
| Some of the drug is dripping out of the needle when you withdraw it from your skin.  | Did you take the needle out of your skin before waiting 5 seconds as directed in Step 15? | If this happens you should inform your doctor. To avoid this problem again, you should always wait 5 seconds after you push the Injection Button before you withdraw the needle from your skin.   |
| The needle is left on the Follistim Pen.   | Have you missed any of the instructions?  | Dispose of the needle in a properly secured container as instructed by your doctor. Change to a new Follistim Cartridge and a new needle.   |
| After your last injection, a remaining volume may be left in the cartridge in addition to the normal quantity of drug dispensed.                   | The cartridge contains extra volume for checking the drug flow.                           | This is not a problem.  |
| You cannot get the cartridge out of the Follistim Pen.   | Is the needle attached?   | Remove the needle from the Follistim Pen and dispose of properly.<br>(Unscrew the Cartridge Holder from the Pen Body and take out the cartridge.)   |
| You are not sure how much drug is left in the cartridge and you do not want to start an injection and then find out that there is not enough drug. | Have you kept good records of your doses?   | In case of any doubt, you should load a new, unused Follistim Cartridge into the Follistim Pen. See “If There is Not Enough Follistim AQ in the Cartridge.”<br><br>To avoid this problem again, you should record your injections. (Women should use a treatment diary.)  |

Important: If you have a question, always mention the Lot number of your Follistim Pen as printed on the Pen Body. If you have a complaint, please do not discard any product or packaging.

For questions on information contained in this leaflet, call 1-866-836-5633.

**www.follistim.com**

### **How Do I Throw Away Used Cartridges and Needles?**

Check with your healthcare provider or pharmacist for instructions about the right way to throw away used cartridges and needles. There may be special local or state laws about how to throw away used syringes and needles.

- **Do not** throw away used cartridges and needles in the household trash and do not recycle them.
- Put used and empty cartridges and needles in a closeable, puncture-resistant container. You may use a sharps container (such as a red bio-hazard container), a hard plastic container with a screw-on cap (such as an empty detergent bottle) or in a metal container with a plastic lid, (such as a coffee can).
- When the container is full, tape around the cap or lid to make sure the cap or lid does not come off.
- When your injection is given by another person, this person must also be careful when removing the cartridge and needle and disposing of the cartridge and needle to prevent accidental needle stick injury and passing infection.

### **How Do I Care for the Follistim Pen?**

1. Clean all exposed surfaces of the Follistim Pen with a clean, damp cloth such as a paper towel. Never wash it in water, detergent or strong medical cleaners.
2. Handle the Follistim Pen carefully to avoid causing damage. You could damage the Follistim Pen by dropping it or handling it roughly.
3. Keep the Follistim Pen away from dust and dirt.
4. Never store the Follistim Pen with a needle attached to it. If you store the Follistim Pen with the needle attached, the drug could leak out and there is risk of contamination.
5. If the Follistim Pen breaks or is damaged, do not try to fix it yourself. Contact your healthcare provider.
6. Do not share your Follistim Pen with another person.

### **How should I store Follistim AQ Cartridge?**

- Store Follistim AQ Cartridge in the refrigerator between 2°-8°C (36°-46°F) until the expiration date.
- Follistim AQ can be stored at or below 25°C (77°F) for 3 months or until the expiration date, whichever comes first. Once the rubber inlay of the Follistim AQ Cartridge has been pierced by a needle, the product may be stored only for a maximum of 28 days at 2°-25°C (36°-77°F).
- Keep Follistim AQ Cartridge away from light.
- Do not freeze.

**Keep Follistim AQ Cartridge, needles, and the disposal container, out of the reach of children.**

Manufactured for: Merck Sharp & Dohme Corp., a subsidiary of  
 **MERCK & CO., INC.**, Whitehouse Station, NJ 08889, USA

Manufactured by: Vetter Pharma-Fertigung GmbH & Co. KG, Ravensburg, Germany

BD, BD Logo and BD Micro-Fine are trademarks of Becton, Dickinson and Company.

For patent information: [www.merck.com/product/patent/home.html](http://www.merck.com/product/patent/home.html)

Copyright © 2004-2019 Merck Sharp & Dohme B.V., a subsidiary of **Merck & Co., Inc.**  
All rights reserved.

Revised: 09/2019

usppi-mk8328-SOI-1909r014